



415190

Dilumat Start



Chemin de l'Orme
69280 MARCY-L'ÉTOILE
FRANCE
RCS LYON 673 620 399
Tel: + 33 (0)4 78 87 20 00
Fax: + 33 (0)4 78 87 20 90

www.biomerieux.com

Document history

Translation of the original document.

VERSION	DATE	MODIFICATIONS
A	27/02/2013	Creation
B	17/09/2013	Modifications on images
C	16/12/2013	Regulatory modifications following the AES and bioMérieux merger

User manual - Pictograms

 Note
 Equipment optimization tip
 Prohibited action notice
 Noteworthy point warning
 Danger or potential risk warning
 Reminder of the pre-requisites for implementation of the next instructions

AES blue line™ is one of bioMérieux's registered trademarks.

The information and graphs contained in this manual are not binding. bioMérieux therefore reserves the right to implement changes to the document without prior notice.

 bioMérieux will not approve any modifications whatsoever to the equipment implemented by the user. bioMérieux shall under no circumstance be held liable for any direct or consequential injury or damage whatsoever suffered by the user or by any third party, as a result of modifications to the equipment.

-  • This equipment is for professional use only.
- Users are required to read all the accompanying documents, including the statutory information, before using the equipment.

Equipment conformity

FCC compliance

This equipment has been declared in compliance with FCC regulations, section 15, applicable to class B digital equipments. These regulations are designed to provide adequate protection against harmful interferences in a residential installation. This equipment generates, uses and may emit radio-electric waves. It can generate interferences that may be harmful to radio-communications if it is not installed or used according to the instructions. Using this equipment in a specific installation may generate harmful interferences, in which case the user may need to correct the interference using one or all of the following methods:

- Redirect or move the receiver aerial
- Move the equipment away from the receiver
- Plug the equipment into a different socket from the one used for the receiver
- Contact the distributor or get help from a radio / TV technician.

Industry conformity Canada (IC)

This class B digital equipment complies with Canadian standard NMB- 03.

EC compliance

This equipment complies with the relevant EC directives and standards listed in the accompanying certificate(s).

Learning about the equipment	8
Features	8
Front view	9
The touch panel	10
Rear view	11
Context-sensitive help	12
Equipment features explained	13
Installing the equipment	14
Installing the tubing	16
Switching the equipment on / off	18
Switching the equipment on	18
Placing the equipment in sleep mode	18
Restarting the equipment	18
Switching off the equipment	18
Recommendations for use	19
Understanding the light indicators	20
Calibrating	23
Making an adjustment	24
Performing a dilution	25
Dosing	26
Dispensing in manual mode	27
Stopping a task under way	28
Recharging diluent while a task is under way	28
Changing the pump while a task is under way	28
Canceling a task under way	28
Priming a tubing	29
Cleaning the equipment	30
Connecting with an administrator account	31
Managing the external devices	32
Enabling-disabling the embedded FTP server	32
Connecting a USB stick	32
Configuring a Wifi stick	33
Configuring the Ethernet	33
Configure a remote FTP server	34
Connecting a physical printer	34
Creating a virtual printer	34
Modifying a label model	35
Importing a label model	35
Configuring the bar code reader	36

Configuring a number pad	36
Configuring the system settings	37
Importing the configuration settings	37
Exporting the configuration settings for the equipment	37
Managing updates	37
Configuring the localization settings	38
Adjusting the energy settings	38
Adjusting the date, time and time zone	38
Configuring the settings for the equipment	39
Print a certificate	39
Restore the calibration and adjustment certificates	39
Configuring the communication settings for the LIMS	39
Managing messages	40
Displaying a message	40
Deleting a message	40
Troubleshooting	41
Typical errors	41
Displaying the details for the equipment	42
Checking sensor state	42
Implementing functional tests	42
Managing the life counters	43
Authorizing remote maintenance	43
Technical data	44
e-pump	45
Technical specifications	46
Recommendations for use	46
Installing an e-pump	47
Glossary	48
Index	49

Learning about the equipment

Features

Diluting the samples is the first step in microbiological analysis. Accuracy in the dilution process therefore determines the success of the whole analysis. The Dilumat Start allows you to optimize this step and to carry it out automatically.

Place the sample in the bag. It will be weighed accurately by the equipment then the exact volume of diluent will be automatically and carefully added to achieve the expected dilution or dosing outcome.

The Dilumat Start gives you the benefit of the latest technological innovations to help make your job much easier.

✓ **New!**

- Get detailed explanations thanks to the Dilumat Smart's integrated context-sensitive help. The features are presented in a way that will guide you throughout your use of the equipment.

✓ **Ergonomics**

- The color touch-screen and the different light indicators make it easier to find what you need.
- The mast light indicator allows you to quickly visualize task results.
- The Dilumat Smart is compact enough to be installed under laminar flow hoods.
- The concept of the dispensing arm allows you to remove the bag to limit the risks of cross contamination.
- The design allows for easy maintenance and an excellent resistance to cleaning products.

✓ **Efficiency and accuracy**

- Diluting a 25 g sample to 1/10 is done in less than 20 seconds.
- The strain gage guarantees 100% accuracy up to 100 g.
- The pump guarantees identical dispensing patterns.


✓ **Capacity and adaptability**

- The equipment can dispense up to 2 kg (sample + diluent).
- You can connect an external pump and manage up to two different diluents on the equipment.

✓ **Advanced features**


- Traceability is guaranteed thanks to data imports and exports via various peripherals (USB stick, ticket or virtual printer, FTP server).
- Traceability of the diluent batches is easily managed thanks to the embedded RFID reader bioMérieux and the bar code reader.

Front view

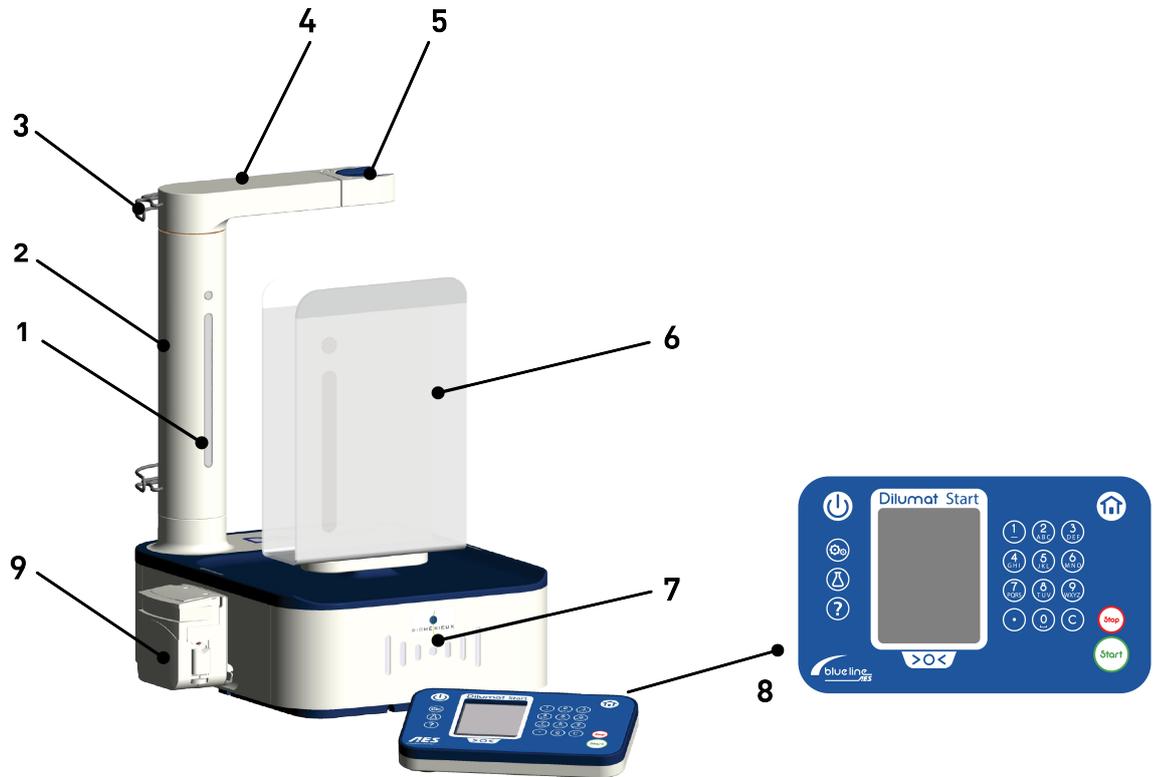


1	Mast light indicator
2	Dispensing mast
3	Hose guide
4	Dispensing arm
5	Dispensing nozzle support
6	Bag support
7	RFID reader with light indicator
8	Touch panel
9	Pump head

The touch panel

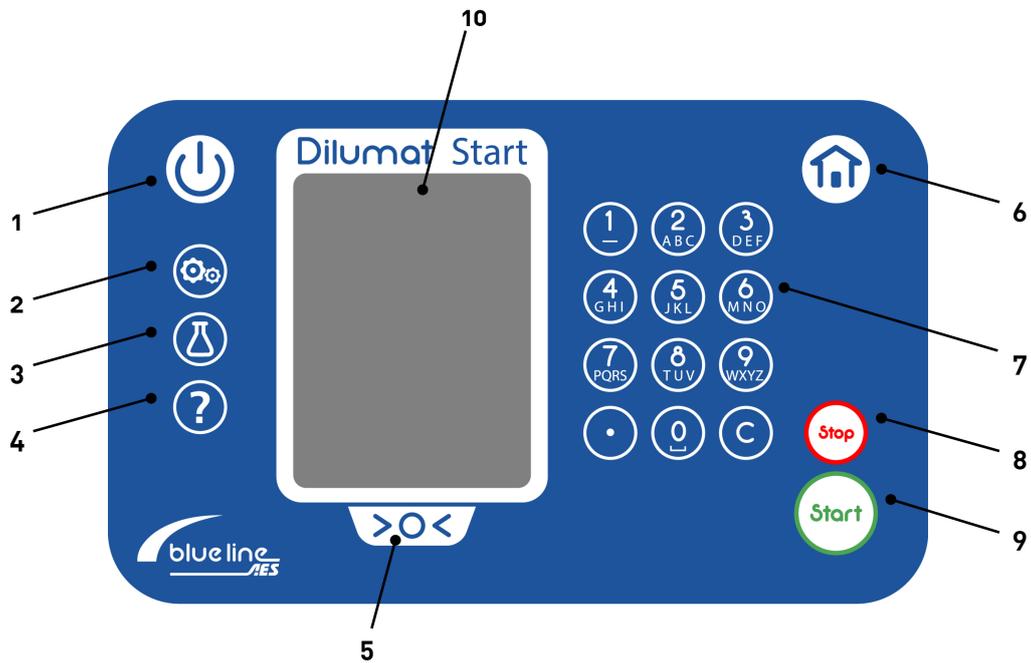

1	On / Off button
2	Configuration button
3	Diluent button
4	Help button
5	Calibration button
6	Home button
7	Number pad
8	Stop button
9	Start button
10	Touch screen

Rear view

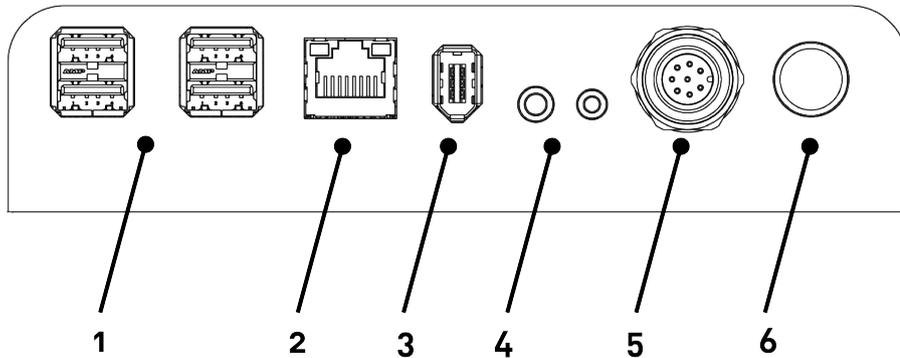

1	USB port x 4
2	Ethernet port
3	Touch panel power supply cable connector
4	Dry input / output contact
5	External pump power supply cable connector
6	Mains power supply cable connector

Context-sensitive help

- Use the context-sensitive help available in the touch panel to find the information you need about specific equipment features.
- Read the user manual to find out how to use the equipment and how to access different functions, then get any additional information you need from the context-sensitive help feature.


Whichever screen you happen to be working on, just press  to display help information on any particular feature.



1. Press.
2. Select.
3. Read.

Just press  to disable the Help feature.

Equipment features explained

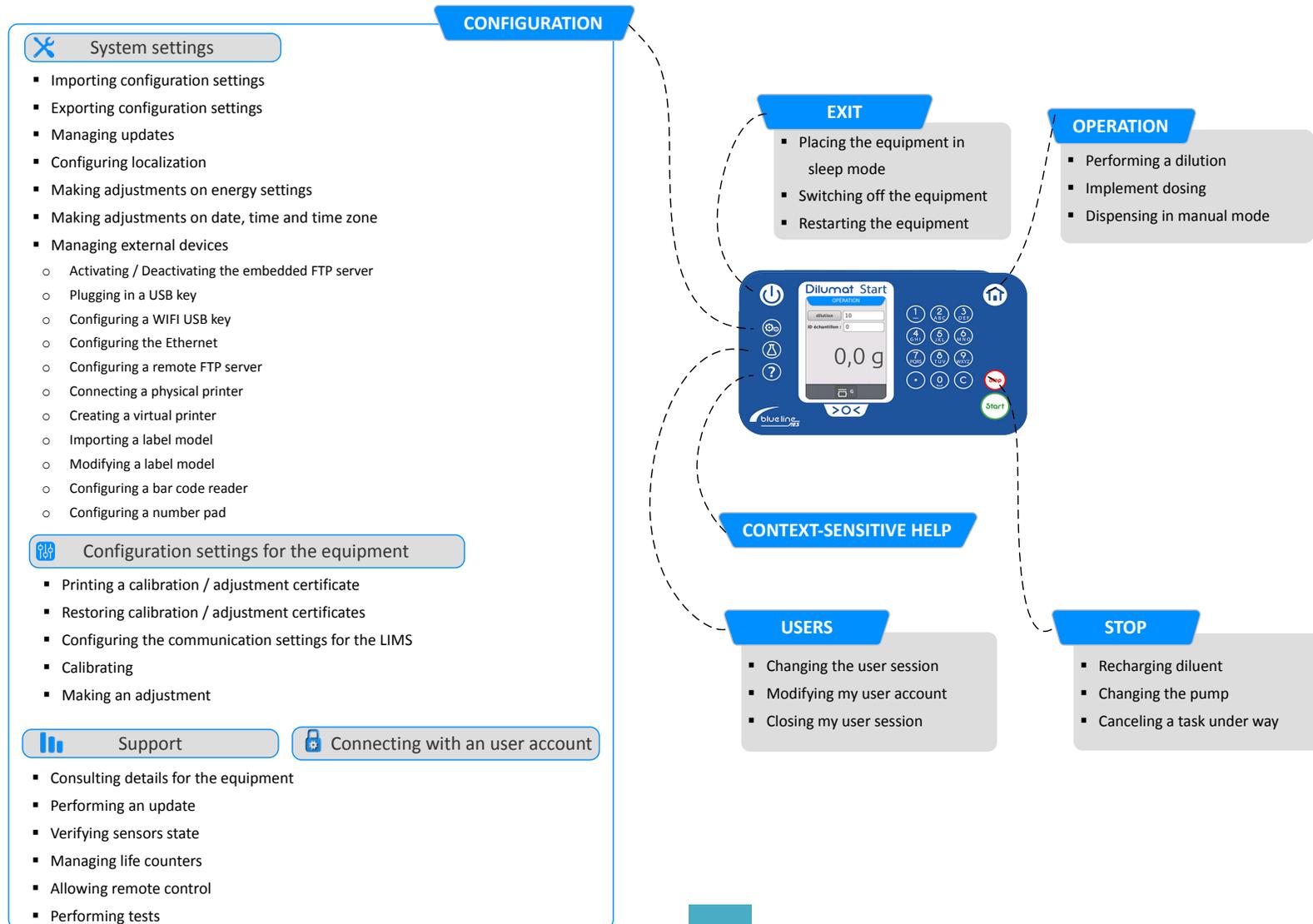

Installing the equipment

Install the equipment in a good working environment for optimal use.











You have unpacked and assembled your equipment according to the assembly procedure.

- Carefully place the equipment on a flat, horizontal, stable surface.
The equipment should be placed on a clean, dry work surface.
- You may run the cable through the grooves under the equipment depending on where you want to place the touch panel.



Do not connect the touch panel before performing this operation.

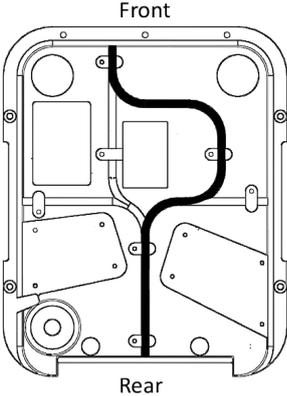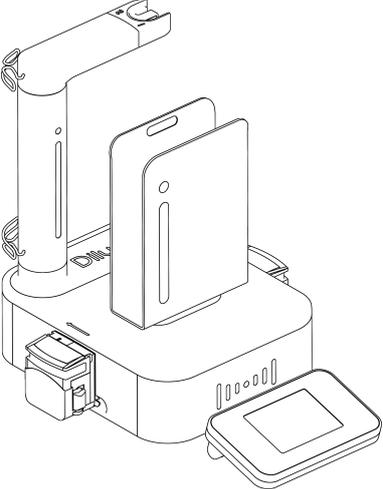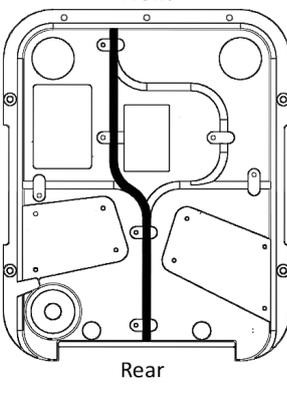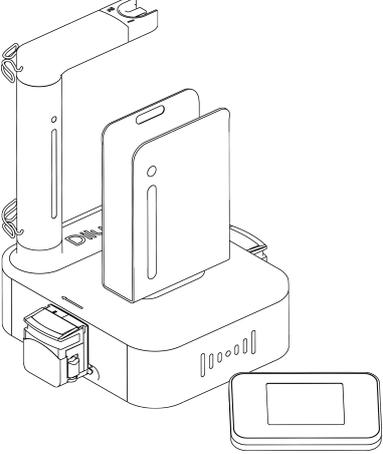
Location of the touch panel	Underneath the equipment	Real preview
In front of the equipment with a short cable length	 <p>Front</p> <p>Rear</p>	
In front of the equipment with a long cable length	 <p>Front</p> <p>Rear</p>	

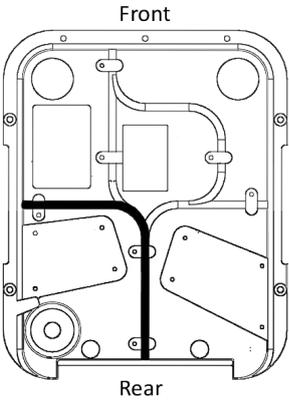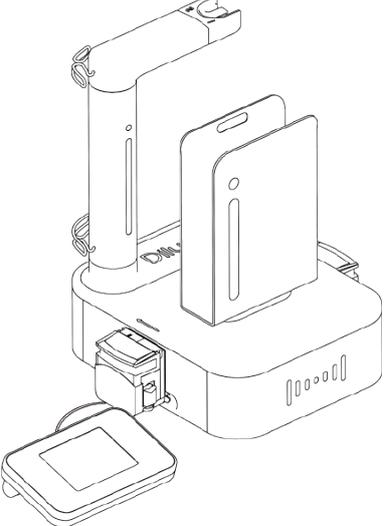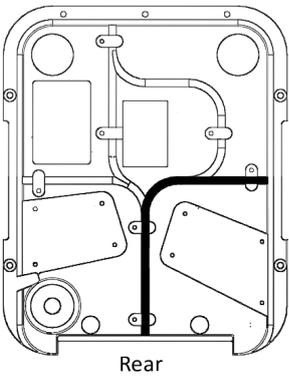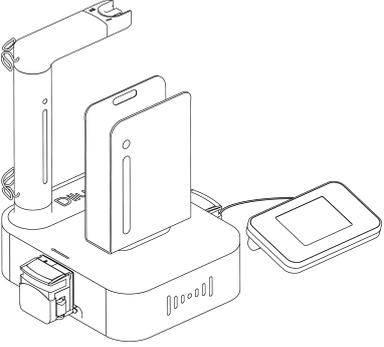Location of the touch panel	Underneath the equipment	Real preview
Left side of the equipment		
Right side of the equipment		


3. Connect the touch panel to the equipment with its power cable.
4. Plug the equipment's power cable into a wall socket close to the device. In the event of a problem, the user must be able to easily unplug the power cable.
The equipment will start up.
5. Choose a country that has localization settings that are similar to the settings in your country. You can then customize the localization settings adapted to your country. Check the "Define the localization settings" section.
The Task screen will appear.

Installing the tubing

1. Lift up the pump head housing.

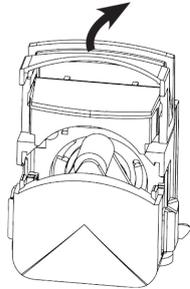


2. Depending on the tubing being used, turn the thumb wheels to adjust the pump head.

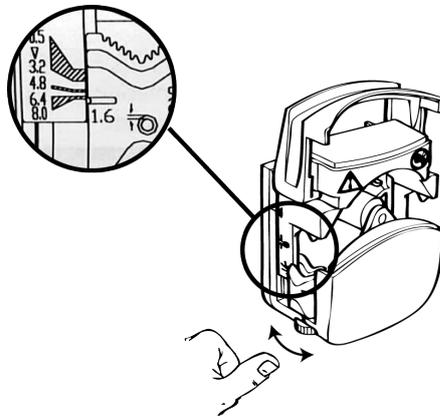

3. Insert the tubing (tip of the opaque tubing) into the pump and make sure that it is placed in the right direction.


The part of the tubing with the dispensing nozzle must be placed in the direction of the arrow.

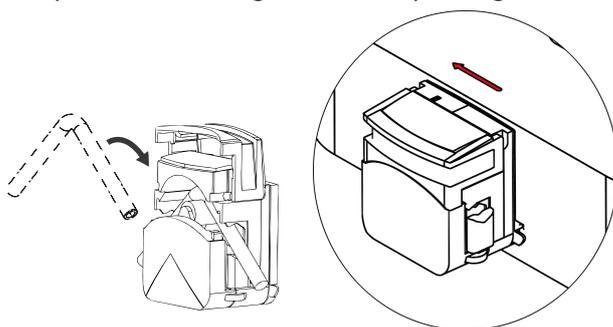


4. Close the housing while making sure you remove your fingers from the pump.



5. Check that the tubing is properly placed.

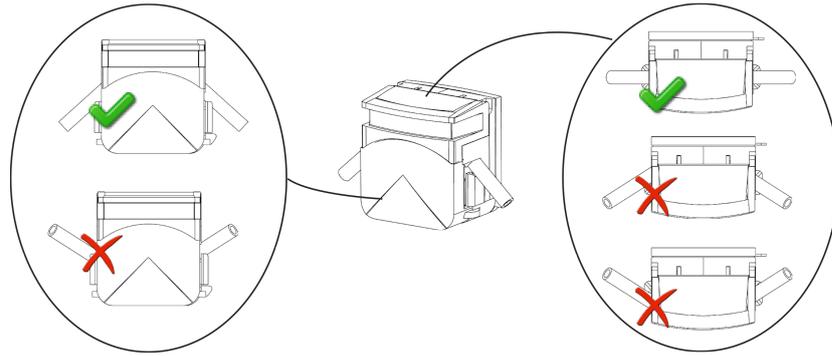


6. Attach the dispensing nozzle to the corresponding nozzle holder while making sure that you do not touch the tip of the nozzle:

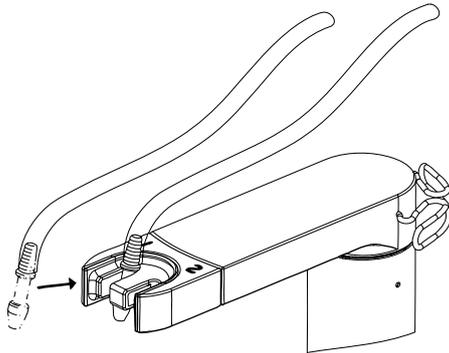

Hose installed on:	Dispensing nozzle placed on:
Left pump	nozzle holder 1
Right pump or external pump	nozzle holder 2

7. Run the tubing, on the dispensing nozzle side, through the tubing guide. Avoid stretching the tubing, doing so might prevent the arm from working properly.

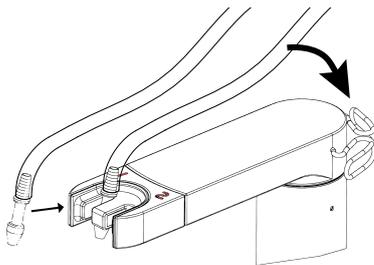

Switching the equipment on / off

Switching the equipment on

 The equipment does not have an on / off switch.

Proceed according to the following situations:

Option	Description
Your equipment is off and plugged in.	Press  on the touch panel.  The equipment will start after a few seconds.
Your equipment is in sleep mode.	Touch the screen.
The equipment is not plugged in.	Plug in the equipment.  The equipment will start up after a few seconds.

Placing the equipment in sleep mode

Place the equipment in sleep mode to save energy.

1. Press .
2. Press .

To exit sleep mode, touch the screen.

Restarting the equipment

1. Press .
2. Press .

Switching off the equipment

1. Press .
2. Press .

Recommendations for use




Be careful to keep your head well away from the rotating arm's reach to avoid harm when moving or handling the equipment.

- Always use the equipment on a horizontal, stable surface.
- Adjust the equipment at least once a year.
- Check that the gap between the bag and the bag support does not exceed 1.5 cm.
- Use the following tubing set dimensions:
 - 4.8 mm tubing set with no sampling nozzle: for sample dilutions or dosing with a net weight of under 100 g
 - 6.4 mm tubing set with no sampling nozzle: for sample dilutions or dosing with a net weight of under 1,000 g

Understanding the light indicators

Different light indicators inform you about the state of the equipment and help you to understand how it operates.

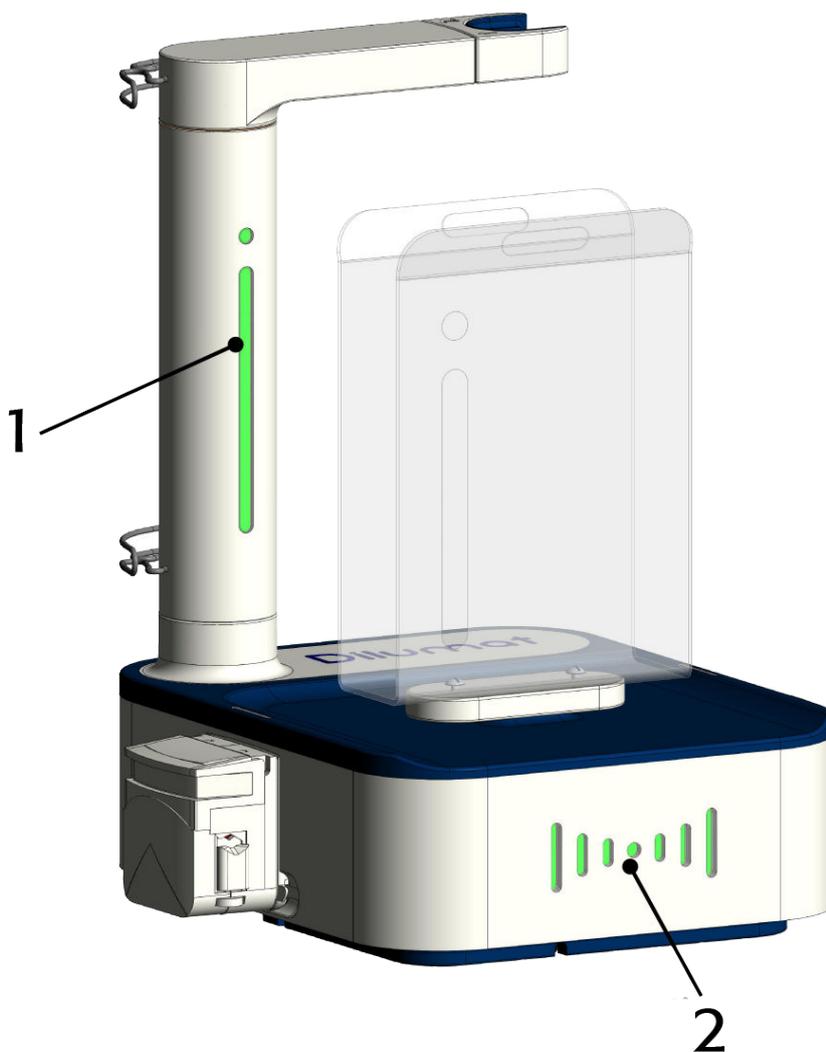

1	Arm light indicator
2	RFID reader light indicator


✓ **Mast light indicator**

Indicator	Color	Explanation
	blue	<ul style="list-style-type: none"> equipment is turned on equipment is ready for an operation
	blue (flashing)	<ul style="list-style-type: none"> task in progress task in manual mode: manual dispensing of the diluent currently in progress
	green	<ul style="list-style-type: none"> task successful
	green (flashing)	<ul style="list-style-type: none"> Dispensing in manual mode: the correct diluent mass has been reached.
	red	<ul style="list-style-type: none"> task failed an error has occurred
	orange (flashing)	<ul style="list-style-type: none"> task postponed
	colorless	<ul style="list-style-type: none"> equipment is off equipment is in standby mode

✓ **RFID indicator**

Indicator	Color	Explanation
	Briefly turns green)	<ul style="list-style-type: none"> label or RFID card scan was successful
	blue	<ul style="list-style-type: none"> RFID reader is in playback mode
	colorless	<ul style="list-style-type: none"> cannot be scanned
	red	<ul style="list-style-type: none"> scan failed

✓ **RFID indicator (with LIMS connection)**

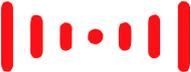Indicator	Color	Explanation
1.  2.  3. 	1. flashing blue 2. turns green 3. going back to blue	<ul style="list-style-type: none"> RFID ticket or RFID card scan was successful
	blue	<ul style="list-style-type: none"> RFID reader is in playback mode
	red	<ul style="list-style-type: none"> scan failed
	colorless	<ul style="list-style-type: none"> cannot be scanned

Calibrating

Use a standard mass to calibrate the equipment and check whether your equipment readings are accurate.

1. Remove all items from the scales
2. Press the  key then the  key.
3. Press .
4. Press  then **Start** to access the Calibration screen.
5. Enter the calibration mass.
6. Place the calibration mass on the scales.
The equipment displays the mass reading.
7. Press **To generate the certificate**.

Making an adjustment

With the help of the calibration weights, you can make an adjustment to the equipment to reduce the gap between the calculated measurement and the calibration weights.

1. Remove all objects from the scales.
2. Press the  key then the  key.
3. Press  then .
4. Press **Start** to go to the Adjustment screen.
5. If necessary, change the calibration weight for the first adjustment.
6. Put the calibration weight on the scale.
The equipment will indicate the weight measurement.
7. Press **Save**.
The results are saved and appear in the table.
8. If necessary, change the calibration weight for the second adjustment.
9. Put the calibration weight on the scales.
The equipment will indicate the calculated weight.
10. Press **Save**.
The results are saved and appear in the table.
11. Press **Generate certificate**.

You have to perform a calibration to confirm the adjustment certificate, check the [Perform a calibration](#) section.

Performing a dilution

Distribute a diluent according to a specific dilution factor and sample mass.

1. Press  to access the Task screen, then select **Dilution** and enter the dilution factor you want to achieve.
2. Press  to select a pump..
3. Press  to allocate a diluent to the pump you have just selected.
4. Place your bag on the bag support.
5. Carry out a tare weight measurement.
6. Place the sample in the bag and complete the **sample ID** field if required.
7. Press  .
Dispensing starts.
8. Remove the bag once the task is complete.

Dispense a certain mass of diluent.

1. Press  to access the Task screen, then select **Dosing** and enter the mass of the dose you need to dispense.
2. Press  to select a pump.
3. Press  to allocate a diluent to the pump you have just selected.
4. Place your bag on the bag support.
5. Carry out a calibration.
6. Place the sample in the bag and complete the **sample ID** field if required.
7. Press  .
Dispensing starts.
8. Remove the bag once the task is complete.

Dispensing in manual mode

Dispensing manually by pouring the diluent into the bag.

1. Press  to access the Task screen, then select **Dilution** or **Dosing** and enter the mass of the dose you need to dispense or the dilution factor you want to achieve.
2. Press  to select the **Manual** dispensing mode.
3. Place your bag on the bag support.
4. Carry out a calibration.
5. Place the sample in the bag and complete the **sample ID** field if required.
6. Press  .
The mast light indicator flashes blue.
7. Pour the diluent into the bag.
8. Press  to confirm.
9. Remove the bag once the task is complete.

Stopping a task under way

You may stop a task under way at any time to top up the diluent or change the pump. You can then resume or cancel the task.

Recharging diluent while a task is under way

You may recharge the diluent while a task is under way, without canceling the task. This only applies if the Diluent option is enabled and if a diluent has been allocated to the pump being used for the task.

A task is under way.

1. Press  .
The Stop screen is displayed.
2. Press  to recharge the diluent. The operation starts automatically.

Changing the pump while a task is under way

You may change the pump while a task is under way, without having to cancel the task. This only applies when at least two pumps have been detected by the equipment.

A task is under way.

1. Press  .
The Stop screen is displayed.
2. Select a pump if required.
3. Press  .
The task starts automatically.

Canceling a task under way

You may cancel a task that is under way.

A task is under way.

1. Press  .
The Stop screen is displayed.
2. Press  then  .

Priming a tubing

Prime the tubing of the pump you are using in case of malfunction.



You have positioned the tubing correctly in the pump but the pump stalls when you start a task.


1. Press  to select the pump to which the tubing has been connected.
2. Press .
3. Press  until the diluent flows out of the pump head.

Cleaning the equipment



- Do not spray cleaner directly on the equipment terminals.
- Do not use aggressive chloride- or acid-based products.

1. Unplug the equipment from the wall socket.
2. Remove the dispensing tubing from the equipment.
3. Remove the pump head, if required. Be careful not to trap your fingers.

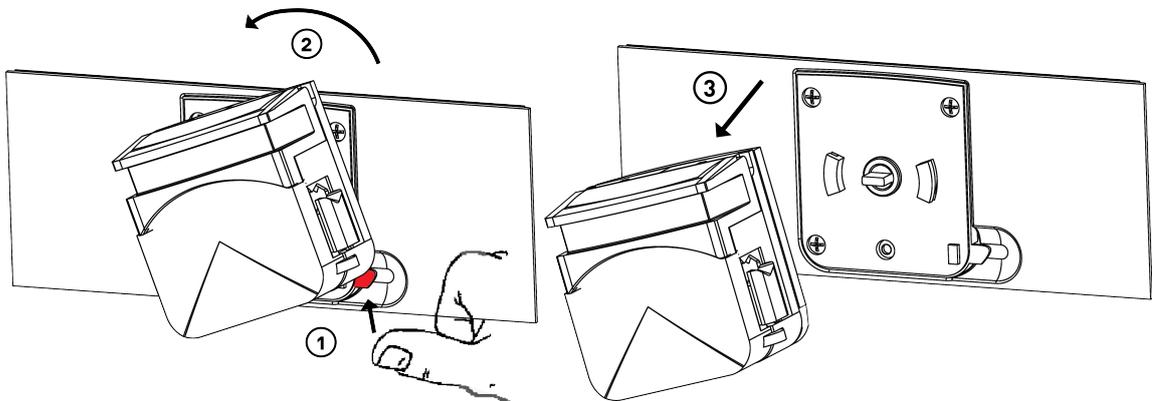

4. Remove the bag support.
5. Clean the equipment, the pump head and the bag support with a cloth and a cleaning fluid or with hot water.
Remember to also clean the inside of the arm.
6. Clean the touch panel with a dry or slightly damp cloth

Connecting with an administrator account

You can connect with an administrator account to access a wider range of features.

1. Press  .
The Configuration screen is displayed.
2. Press  .
3. Enter the password 35170 to connect with an Administrator account.
You will remain connected as an Administrator until you leave the Configuration screen.

Managing the external devices

Access more features (data management, updates, etc.) by adding physical or virtual external devices (USB stick, Wifi stick, Ethernet, printer, virtual printer, number pad, bar code reader, embedded FTP, remote FTP).

Enabling-disabling the embedded FTP server

You may enable or disable the embedded FTP server to store your data.

1. Press the  key then the  key.
2. Press .
3. Press .
4. Select the embedded FTP in the list of external devices.
5. Press the on / off button to enable or disable the embedded FTP server.

 You may connect to the embedded FTP server via an anonymous session by using the 35170 password.

Connecting a USB stick

You may use a USB stick to store your data.

 You can connect your USB stick from any screen.

1. Press the  key then the  key.
2. Press .
3. Press .
4. Plug your USB stick into the rear panel on the equipment.
The USB stick is automatically detected and displayed in the list of external devices.

 Disable the USB stick before removing it to avoid data loss.

Configuring a Wifi stick

Use a Wifi stick to access a network and to manage your data.



Disable the Ethernet before enabling the Wifi stick.

1. Press the  key then the  key.
2. Press .
3. Press .
4. Plug your Wifi stick into the rear panel on the equipment.
The Wifi stick is automatically detected and displayed in the list of external devices.
5. Press the on / off button to enable the Wifi stick.
6. Press  to configure the Wifi stick.

Configuring the Ethernet

You may use the Ethernet to access a network and to manage your data.



Disable the Wifi stick before enabling the Ethernet.

1. Press the  key then the  key.
2. Press .
3. Press .
4. Plug your Ethernet cable into the rear panel on the equipment.
The Ethernet is automatically detected and displayed in the list of external devices.
5. Press the on / off button to enable the Ethernet.
The equipment connects to the network.
6. Press  to display the IP settings.
By default, the equipment is configured in DHCP mode.

Configure a remote FTP server

You may configure a remote FTP server to store your data.

1. Press the  key then the  key.
2. Press .
3. Press .
4. Press .
5. Select **remote FTP server** then press .
6. Complete the fields, then press  to confirm.
The remote FTP server is included in the list of external devices.

Connecting a physical printer

You may use a physical printer to print your labels.

1. Press the  key then the  key.
 2. Press .
 3. Press .
 4. Plug your printer into the rear panel on the equipment.
The printer is automatically detected and displayed in the list of external devices.
-  The virtual printer is configured by default.
5. Press  to modify its configuration.

Creating a virtual printer

You may use a virtual printer to store your data on a USB stick, an embedded FTP server or a remote FTP server.

1. Press the  key then the  key.
 2. Press .
 3. Press .
 4. Press .
 5. Select **virtual printer** then press .
 6. Complete the fields, then press  to confirm.
The virtual printer is displayed in the list of external devices.
-  The virtual printer is configured by default.
7. Press  to modify its configuration.

Modifying a label model

You may modify the default label models for the equipment.



Your printer is on or you have set up a virtual printer.

1. Press the  key then the  key.
2. Press .
3. Press .
4. Select a printer and press .
5. Select the type of label you want to modify.
6. Choose another label model.
7. Press  to confirm.

Importing a label model

You may customize the label models by importing your own models.



- Your printer is on or you have set up a virtual printer.
- You have stored your label model file on an external storage device (USB stick, embedded FTP server or remote FTP server).

1. Press the  key then the  key.
2. Press .
3. Press .
4. Select a printer and press .
5. Select the type of label you want to customize.
6. Press .
7. Select and import your label model file from the external storage device.

Configuring the bar code reader

You may use a bar code reader to scan your labels.

1. Press the  key then the  key.
2. Press .
3. Press .
4. Plug your bar code reader into the rear panel on the equipment.
The bar code reader is automatically detected and is included in the list of external devices.
5. Press  to configure the bar code reader.

 A bar code reader is configured like a keyboard, according to the country's language. For example in France (fr), keyboards are configured in AZERTY and in QWERTY in English speaking countries. It is the same for bar code readers. Be sure to select the country code corresponding to your bar code reader localization settings. Otherwise, the reader will not be able to read correctly your bar codes. Run a test on a bar code on which you can see the content in order to check that your reader has been correctly configured.

6. Select your country's language
7. Press  to confirm.

Configuring a number pad

You may use a number pad.

1. Press the  key then the  key.
2. Press .
3. Press .
4. Plug your number pad into the rear panel on the equipment.
The number pad is automatically detected and is included in the list of external devices.
5. Press  to configure the number pad.
A number pad is configured like a keyboard, according to the country's language.
6. Press  to confirm.

Configuring the system settings

Importing the configuration settings

You may import your configuration settings for the equipment from your external devices.

1. Press the  key then the  key.
2. Press .
3. Press  to access the Import/Export screen.
4. Press .
5. Select the configuration file you want to import.
6. Press  to confirm.
The equipment will restart.

Exporting the configuration settings for the equipment

You may export your configuration settings for the equipment via your external devices.

1. Press the  key then the  key.
2. Press .
3. Press  to access the Import / Export screen.
4. Press .
5. Select where you want to store the configuration file.
6. Press  to confirm.

Managing updates

You may enable and disable the automatic firmware installation and search feature and update the firmware.

1. Press the  key then the  key.
2. Press  then  to access the Update screen.
3. Configure the settings.
4. Press  to confirm.

Configuring the localization settings

You may choose your localization settings by choosing a country.

1. Press the  key then the  key.
2. Press .
3. Press  to access the Localization screen.
4. Select a country.
5. Press  to confirm.

Adjusting the energy settings

You may adjust the energy settings for greater readability and configure the sleep feature to save energy.

1. Press the  key then the  key.
2. Press .
3. Press  to access the Energy screen.
4. Configure the settings.
5. Press  to confirm.

Adjusting the date, time and time zone

You may adjust the date, time and time zone You may also enable / disable the automatic date / time and time zone synchronization feature.

1. Press the  key then the  key.
2. Press .
3. Press  to access the Date & Time screen.
4. Complete the fields
5. Press  to confirm.

Configuring the settings for the equipment

Print a certificate

You may print the calibration and adjustment certificates using a physical or a virtual printer. The labels are printed as .png files and can be stored on a USB stick, an FTP server, etc.



Check that you have configured all the external devices required for printing. See the [Managing your external devices](#) section.

1. Press the  key then the  key.
2. Press .
3. Press  to access the Certificates screen.
4. Select the benchmarking or adjustment tab according to the certificate you want to print.
5. Press  to print the certificate.
The certificate will be printed as a label or saved as a .png file on a USB stick or FTP server.

Restore the calibration and adjustment certificates

You may restore the default values of the calibration and adjustment certificates.

1. Press the  key then the  key.
2. Press .
3. Press  to access the Certificates screen.
4. Press  to restore the default values of the calibration and adjustment certificates.

Configuring the communication settings for the LIMS

You may configure the settings to allow communication with your LIMS.



Contact bioMérieux customer service to configure the LIMS.

1. Press the  key then the  key.
2. Press .
3. Press .
4. Configure the settings.
5. Press  to confirm.

Displaying a message

A flashing envelope indicates you have a new message, no matter which screen you are working on.

1. Press .
The Message screen is displayed.
2. Select the message you want to read, then press .
 - The message is marked as read.
 - It remains in your inbox.
 - The envelope stops flashing.

Deleting a message

Delete your messages from the message list once you have read them.

1. Press , no matter which screen you are working on.
The Message screen is displayed.
2. Select the message you want to delete, then press .

Typical errors

Please check the following trouble-shooting list before calling the after-sales service.

Problem outline	Possible causes	Problem resolution
Pump will not start	<ul style="list-style-type: none"> Badly placed tubing Hose not primed 	<ul style="list-style-type: none"> Place tubing correctly in pump Prime tubing.
Uneven dispensing.	<ul style="list-style-type: none"> Tube adjustment is not in right position for the diameter of the tubing currently in use. 	<ul style="list-style-type: none"> Adjust the tube so that the tubing is not constricted.
The touch panel is off but the mast light indicator is on.	<ul style="list-style-type: none"> Faulty touch panel connection on the equipment side. Faulty touch panel connection on the touch panel side. 	<ul style="list-style-type: none"> Check that the touch panel cable is connected to the equipment. <ol style="list-style-type: none"> Disconnect the touch panel from the equipment. Unscrew the touch panel connector housing. Re-connect the cable correctly. Screw on the touch panel connector housing. Re-connect the touch panel to the equipment.
The mast light indicator is off but the touch panel is on.	<ul style="list-style-type: none"> Faulty mast connection. 	<ul style="list-style-type: none"> Check that the mast cable is connected to the terminal in the base of the equipment.
Only half the touch panel screen is on.	<ul style="list-style-type: none"> Faulty connection 	<ul style="list-style-type: none"> Disconnect then re-connect the equipment.
Dispensing will not start.	<ul style="list-style-type: none"> Dispensing in manual mode. 	<ul style="list-style-type: none"> Select the required dispensing mode. In manual dispensing mode, you must introduce the diluent manually, without using the pump.
Weight remains the same.	<ul style="list-style-type: none"> The equipment has reached its maximum weighing capacity. 	<ul style="list-style-type: none"> Check that total weight (sample + diluent) does not exceed the equipment's weighing capacity.
Equipment is not dispensing diluent in the middle of the bag support.	Faulty tubing nozzle positioning	<ul style="list-style-type: none"> Check that all tubing nozzles are correctly positioned on the end connectors.
Dispensing starts but no liquid is being dispensed.	<ul style="list-style-type: none"> Hose is constricted. Hose is positioned wrong way round. 	<ul style="list-style-type: none"> Make sure the tubing is not constricted or bent. Make sure the tubing are correctly positioned on the left-hand and right-

Problem outline	Possible causes	Problem resolution
		hand pumps by checking the arrows on the equipment.
Bar code reading malfunction.	<ul style="list-style-type: none"> • Incorrect bar code reader configuration. 	<ul style="list-style-type: none"> • Configure the bar code reader when you connect the equipment.
The external pump is not detected on the interface and none of the light indicators are on.	<ul style="list-style-type: none"> • Faulty connection 	<ul style="list-style-type: none"> • Check that the external pump is connected to the equipment.
Pedal does not work.	<ul style="list-style-type: none"> • Faulty connection 	<ul style="list-style-type: none"> • Check that the pedal is connected to the equipment then check state of the pedal entry sensor. Read the Checking sensor state section.
USB ports do not work.	USB devices overload	<ol style="list-style-type: none"> 1. Turn off the equipment. 2. Unplug the equipment. 3. Plug back the equipment.

Displaying the details for the equipment

Change the equipment name, display the equipment specifications and program release.

1. Press the  key then the  key.
2. Press .
3. Press .
4. If required, change the equipment name and display the available information.

Checking sensor state

Check that equipment sensors are operating correctly.

1. Press the  key then the  key.
2. Press  then .
3. Follow context-sensitive help to check state of sensors.
When state is **on**, the sensor is operational.

Implementing functional tests

You may implement simple tests (light indicators, pumps, sound, etc.) and sequential tests on the equipment to check for malfunction.

1. Press the  key then the  key.
2. Press  then .

Managing the life counters

You may display and zero the life counters for the equipment and send anonymous reports to optimize them.

1. Press the  key then the  key.
2. Press  then .
3. You can either:
 - Display the life counters
 - Enable / disable the sending of regular anonymous reports to bioMérieux.
4. Press  to confirm.

Authorizing remote maintenance

Authorize remote maintenance in the event of equipment failure to allow the After-sales service to take control of your equipment remotely to deal with the problem.



To authorize remote maintenance, you need to be in contact with the After-sales service.

1. Press the  key then the  key.
2. Press  then .
3. Press .
4. Follow the After-sales service instructions.

Technical data

Environmental conditions	Altitude	up to 2,000 m	
	Temperature	from 0 to 45°C max.	
	Relative humidity (RH)	from 10 to 80%	
	Mains supply voltage	100-240 CAV, 50/60 Hz, 2 A	
	Machine supply voltage	24 V DC, 3,75 A	
	Degree of pollution	2	
Atmospheric pressure		700 hPa to 1,100 hPa	
Power rating		90 W	
Installation Category		Type II according to CEI664 directive	
Noise level		67 db	
Equipment dimensions (W x H x D)		Equipment alone (mm) <ul style="list-style-type: none"> Start: 324 x 437 x 333 Expert: 377 x 437 x 333 Expert Evo: 377 x 527 x 336 	Packaged equipment (mm) <ul style="list-style-type: none"> Start: 580 x 330 x 400 Expert: 580 x 330 x 400 Expert Evo: 580 x 346 x 400
Weight of the equipment		Equipment alone (kg) <ul style="list-style-type: none"> Start: 11.6 Expert: 13 Expert Evo: 14 	Packaged equipment (kg) <ul style="list-style-type: none"> Start: 15.5 Expert: 17 Expert Evo: 18
Dilution factor		<ul style="list-style-type: none"> Start: From 1/2 to 1/100th Expert and Expert Evo: From 1/2 to 1/1000th 	
Dispensing speed		Less than 10 seconds for a 1/10 dilution of a 25 g sample.	
Accuracy		<ul style="list-style-type: none"> Weight: less than 1% starting from 5 g Dilution / dosing: less than 5% on the final weight starting from 10 g. 	
TFT resistive touchscreen		Size <ul style="list-style-type: none"> Start: 8.89 cm (3.5 inches) Expert and Expert Evo: 12.70 cm (5 inches) 	Resolution (px) <ul style="list-style-type: none"> Start: 240 x 320 Expert and Expert Evo: 640 x 480
Keypad		<ul style="list-style-type: none"> Start: 20 physical keys Expert and Expert Evo: 10 physical keys 	
Weight		<ul style="list-style-type: none"> Start: from 3 g to 2 kg Expert: from 0.3 g to 3 kg Expert Evo: from 3 g to 7 kg 	
External pump	Number of external pumps allowed	<ul style="list-style-type: none"> Start: 1 depending on equipment configuration Expert & Expert Evo: up to 4 	

e-pump



✓ **New**

You can now connect external pumps to your compatible AES blue line™ equipment.

✓ **Capacity**

You can connect up to four external pumps according to the way your equipment is configured.

✓ **User-friendly design**

Different colored light indicators make the system easy to use. Thanks to its compact, robust design, the system allows you to stack the pumps on top of each other.

Technical specifications

Environmental conditions for external pump use and technical specifications.

Environmental conditions

Name	Values
Altitude	up to 2,000 m
Temperature	0°C to 45°C max.
Relative humidity (RH)	10% to 80%
Can be interconnected with:	compatible AES blue line™ equipment
Pollution level	2

Name	Values
Atmospheric pressure	700 hPa to 1100 hPa
Equipment category	Type II ECI664 directive compliant
Noise level	63 db
Dispensing speed	up to 600 rpm

Dimensions and weight

Name	Values
Equipment dimensions (unpacked) (W x H x D)	120 x 115 x 195 mm
Equipment dimensions (packed) (W x H x D)	285 x 340 x 127 mm
Power cable length	400 mm
Power cable extension length	700 mm max.
Equipment weight (unpacked)	3.4 kg
Equipment weight (packed)	3.7 kg

Recommendations for use

- Only use the pump with a single pump head only.
- Total length of power supply cables for external pumps must not exceed 3 m.

Installing an e-pump



1. Place the pump on a flat, stable surface next to the equipment you want to connect it to.
The pump must be placed on a clean, dry work surface.
2. Switch the equipment off and unplug the mains power supply for the equipment you want to connect the external pump to.
3. Connect the power supply cable for the external pump to the equipment.
Connect the second external pump to the pump already connected to the equipment.
4. Plug in the equipment.
 - A blue light indicator will display the number of the external pump according to the order in which they are connected.
 - The green light indicator shows that the external pump has been selected for a task via the equipment.



Adjustment

Adjustment consists in setting the equipment to reduce the difference between the measurement recorded and the reference mass.

Calibration

Calibration has to deal with comparing the measurement recorded by the equipment with a reference mass.

Dilution

Dispensing a diluent according to the dilution factor and sample mass.

Dosing

Dispensing a fixed volume of diluent.

LIMS

LIMS (Laboratory Information Management System), is an integrated management software package used in laboratories.

MPD (Maximum Permissible Difference)

In metrology, this has to deal with the maximum permissible error for a measuring instrument. The metrological analysis is declared non-compliant if the measurement is outside of the zone defined by the MPD.

RFID

RFID (Radio Frequency Identification), is a method for storing and remotely retrieving data by using devices called "RFID tags" or "RFID transponders".

Test weight

Weight used for adjusting and checking.

A**Adjust**

date	38
energy saver	38
time	38

Adjustment

implement	24
-----------------	----

Administrator account

Connecting	31
------------------	----

Authorize

remote maintenance	43
--------------------------	----

B**Bar code reader**

configure	36
-----------------	----

C**Calibration**

implement	23
-----------------	----

Cancel

task	28
------------	----

Certificate

print	39
restore	39

Change

pump	28
------------	----

Check

state of sensors	42
------------------------	----

Clean

equipment	30
-----------------	----

Configure

bar code reader	36
Ethernet	33
LIMS communication settings	39
localization settings	38
number pad	36
remote FTP server	34
wifi stick	33

Connect

printer	34
USB stick	32

Connecting

administrator account	31
-----------------------------	----

Create

virtual printer	34
-----------------------	----

D**Date**

adjust	38
--------------	----

Define

time zone	38
-----------------	----

Delete

message	40
---------------	----

Details for the equipment

display	42
---------------	----

Diluent

recharge	28
----------------	----

Dilution

perform	25
---------------	----

Disable

embedded FTP server	32
---------------------------	----

Dispensing

manual mode	27
-------------------	----

Dispensing position

test	42
------------	----

Display

details for the equipment	42
equipment recommendations for use	19
equipment technical specifications	44
external pump technical specifications	46
external pump use recommendations	46
help	12
life counters	43
message	40
trouble-shooting advice	41

Dosing

implement	26
-----------------	----

Dry contact

test	42
------------	----

E**Embedded FTP server**

disable	32
enable	32

Enable

embedded FTP server	32
---------------------------	----

Energy saver	
adjust	38
Equipment	
clean	30
equipment in sleep mode	18
install	14
restart	18
switch off	18
switch on	18
Equipment configuration	
export	37
import	37
Equipment in sleep mode	
equipment	18
Equipment recommendations for use	
display	19
Equipment technical specifications	
display	44
Ethernet	
configure	33
Export	
equipment configuration	37
External pump	
install	47
External pump technical specifications	
display	46
External pump use recommendations	
display	46
F	
Firmware	
update	37
H	
Help	
display	12
Hose	
prime	29
I	
Implement	
adjustment	24
calibration	23
dosing	26
Import	
equipment configuration	37
label model	35
install	
external pump	47

Install	
equipment	14
tubing	16

L

Label	
print	39
Label model	
import	35
modify	35
Life counters	
display	43
light indicators	
test	42
LIMS communication settings	
configure	39
Localization settings	
configure	38

M

Manage	
USB external device	32
Manual mode	
dispensing	27
Message	
delete	40
display	40
Modify	
label model	35

N

Number pad	
configure	36

P

Perform	
dilution	25
Prime	
tubing	29
Print	
certificate	39
label	39
Printer	
connect	34
Pump	
change	28
test	42

R

Recharge	
diluent	28
refuse	
remote maintenance	43
Remote FTP server	
configure	34
Remote maintenance	
authorize	43
refuse	43
Restart	
equipment	18
Restore	
certificate	39

S

Sound	
test	42
State of sensors	
check	42
Switch off	
equipment	18
Switch on	
equipment	18

T

Task	
cancel	28
Test	
dispensing position	42
dry contact	42
light indicators	42
pump	42
sound	42
Time	
adjust	38
Time zone	
define	38
Trouble-shooting advice	
display	41
tubing	
install	16

U

Update	
firmware	37
USB external device	
manage	32
USB stick	
connect	32

V

Virtual printer	
create	34

W

Wifi stick	
configure	33



bioMérieux SA

Chemin de l'Orme
69280 MARCY-L'ÉTOILE
FRANCE
RCS LYON 673 620 399
Tel: + 33 (0)4 78 87 20 00
Fax: + 33 (0)4 78 87 20 90

www.biomerieux.com

CE



BIOMÉRIEUX